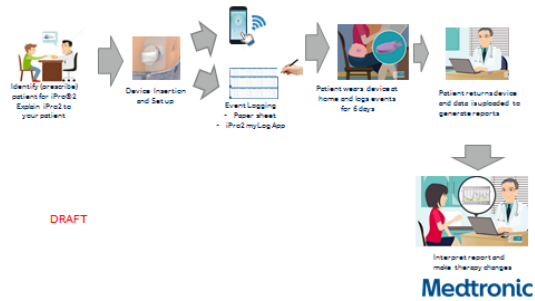

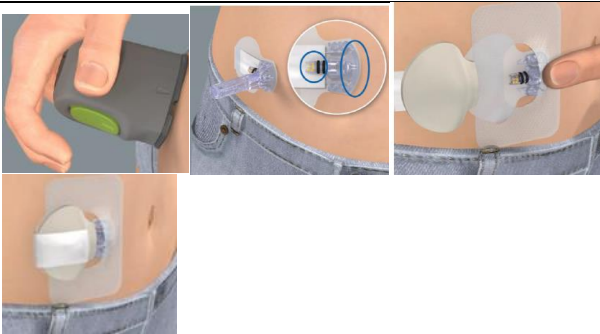



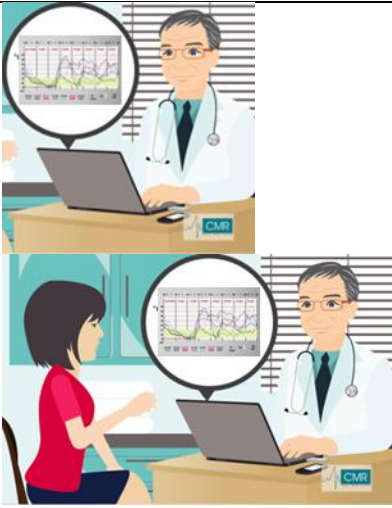
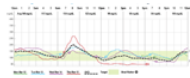


Motion	On-Screen Text	Voice-Over
Show computer with carelink page		Once you have registered your clinic in CareLink® iPro software
Show recorder in docking station		And set up the iPro2 recorder,
<p>Full Steps shown</p> <p>iPro®2 CGM Evaluation Steps</p>  <p>DRAFT</p>	(see slide left for wording/graphics and positioning will change)	You are now ready to start learning about the steps involved in using iPro2 with patients.
 <p>(make sure the patient is clearly an adult)</p>	<p>Identify patient for iPro2 CGM</p> <p>Explain iPro2 CGM to your patient</p>	First, determine or identify which patients will benefit from iPro2 CGM and explain the evaluation process.
 <p>showserter on body, sensor inserted, recorder connected (use iPro recorder with “i”) and then patient with complete sensor and recorder inserted.</p>	<p>Insert sensor</p> <p>Attach recorder</p>	...to place an iPro2 simply insert a sensor into your patient, and then attach the iPro2 recorder.
 <p>Patient with iPro2 device wearing it actively. [ideally we would like to show a patient swimming, as it will provide us with two purposes (watertight and exercise).</p>	<p>Live life as normal</p> <p>Event Logging</p> <ul style="list-style-type: none"> - Paper sheet - iPro2 myLog App 	The patient will be evaluated for up to 6 days and during this time they will record their blood sugars, meals, treatment and exercise.

 <p>Patient using a smartphone vs. patient writing down with a pen on a piece of paper (show them side by side)</p>		<p>This can be done either on a paper logsheet, or using the iPro2 myLog app</p>
<p>Show devices (Sensor and recorder) if possible a patient returning it.</p>	<p>Patient returns device and data is uploaded to generate reports</p>	<p>After the evaluation is over, your patient must return the device to the office along with their paper logsheet if applicable.</p>
<p>Motion of uploading a recorder into the dock</p>		<p>Clean and upload the recorder into the software</p>
 <p>clinician using the computer</p>		<p>If the patient used the app, the logged events will automatically transfer into CareLink iPro software and will be displayed in the report. If not, then transcribe the events from their logsheet.</p>
	<p>Interpret report</p> <p>Make therapy changes</p>	<p>Review the reports generated, and make the desired changes to your patient's therapy regimen.</p>



iPro2 CGM
UNCOVER MORE IN LESS TIME

Medtronic

(see slide left for
wording/graphics and
positioning will change)

Help your patients achieve better
glucose control. Use iPro2 CGM
to uncover more in less time.